


OSKA®

Series3 Auto

Auto-weigh alternating air mattress

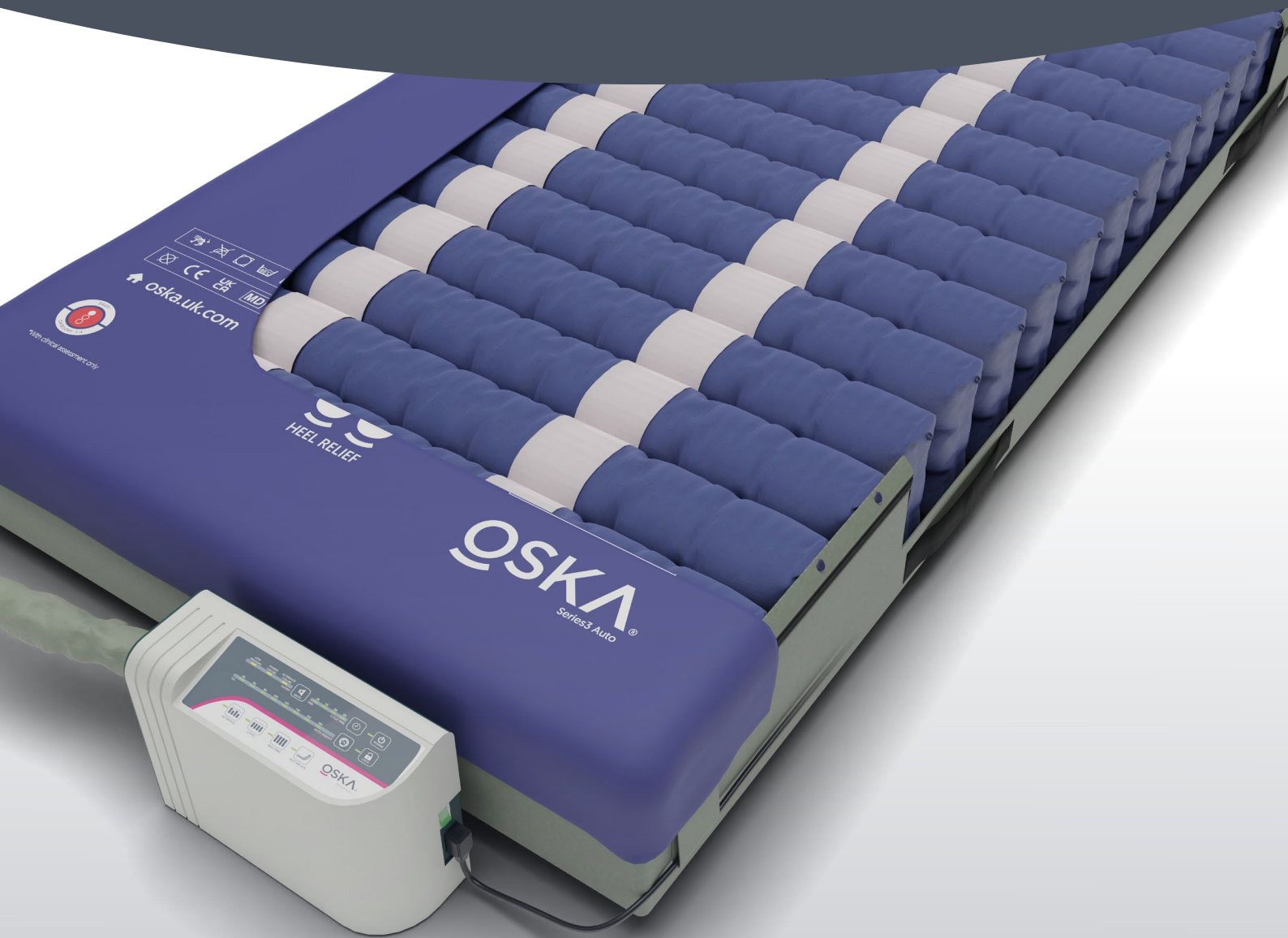

Alternating air, full replacement pressure care mattress.




Auto weigh



Power cut protection



Rapid repair feature

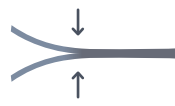

Welded seams



Zoned Offloading

Series3 Auto



Auto-weigh, active therapy alternating air mattress.

The Series3 Auto is a full air mattress replacement system suitable for use across all care settings.

It delivers active therapy to patients and is effective for the prevention of pressure ulcers in patients up to and including those at very high risk of pressure related tissue injury. It can also be used for managing patients with any category of pressure ulcer.

Effective pressure ulcer prevention & management

The alternating 1-in-2 cell cycle achieves regular periods of pressure redistribution and tissue offloading across the patient/support surface interface.


The auto-weigh pump offers a range of therapy settings and automatically adjusts to the patients' weight to ensure optimum air cell pressure at all times.

Bespoke zero pressure zones can also be created by using the unique 'quick connectors' to deflate specific air cells.



Infection prevention

The combination of waterproof cover and welded seams support infection prevention by reducing the risk of cover strike-through and fluid ingress into the mattress.

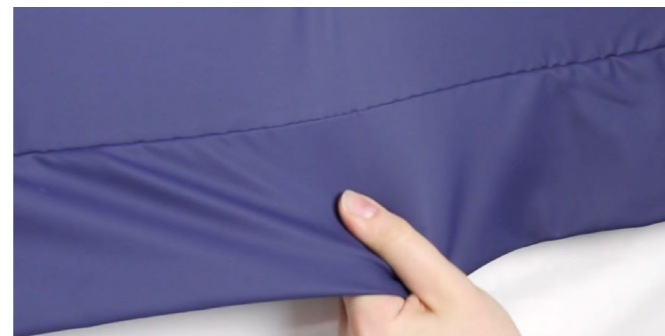


Enhanced patient safety & operational efficiency

The permanently inflated lower portion of each air cell reduces the risk of patients bottoming out during a power cut.

'Quick connectors' facilitate isolation and removal of damaged air cells without interrupting pressure area care therapy. This eliminates 'down-time' due to air cell damage and enhances efficiency by enabling individual air cell replacement in under a minute.



Series3 Auto

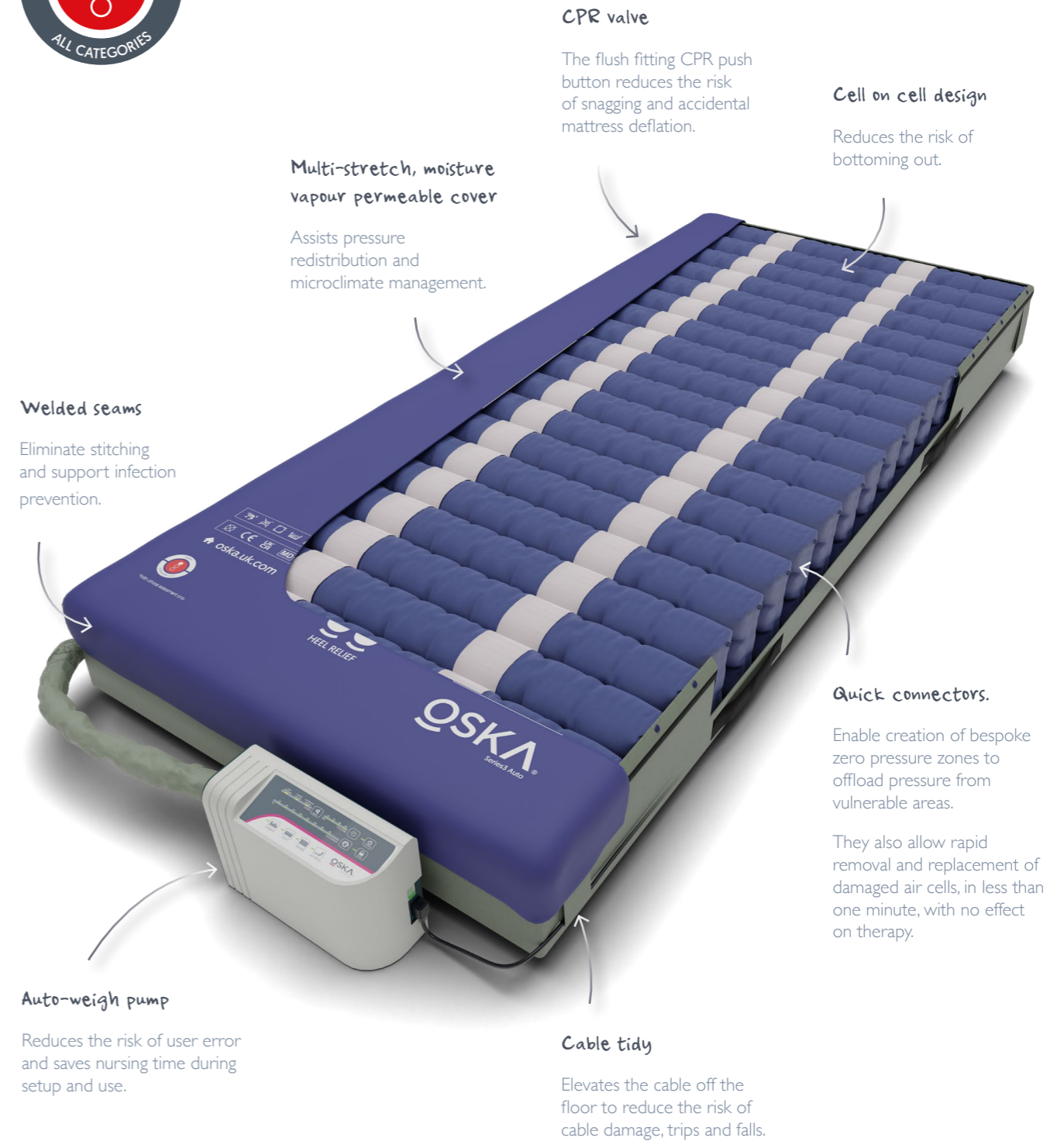

CPR valve

The flush fitting CPR push button reduces the risk of snagging and accidental mattress deflation.

Cell on cell design

Reduces the risk of bottoming out.

Multi-stretch, moisture vapour permeable cover

Assists pressure redistribution and microclimate management.

Welded seams

Eliminate stitching and support infection prevention.

Quick connectors.

Enable creation of bespoke zero pressure zones to offload pressure from vulnerable areas.

They also allow rapid removal and replacement of damaged air cells, in less than one minute, with no effect on therapy.

Auto-weigh pump

Reduces the risk of user error and saves nursing time during setup and use.

Cable tidy

Elevates the cable off the floor to reduce the risk of cable damage, trips and falls.

Also available as:

Series3 Bariatric

Active therapy support surface.



AVAILABLE FOR RENTAL

CPR Valve

The flush fitting CPR push button reduces the risk of snagging and accidental mattress deflation.

Welded seams

Eliminate stitching and support infection prevention.

Multi-stretch, moisture vapour permeable cover

Assists pressure redistribution and microclimate management.

Cable tidy

Elevates the cable off the floor to reduce the risk of cable damage, trips and falls.

Cell-on-cell design

Reduces the risk of bottoming out.

Quick connectors

Quick connectors enable creation of bespoke zero pressure zones to offload pressure from vulnerable areas.

They also allow rapid removal and replacement of damaged air cells in less than a minute with no effect on therapy.

Digital pump*

Simple to set up and use, with manually selected weight settings.

*N.B. Series3 Bariatric does not have an autoweigh pump.



Specification

Prevention	Very High Risk*
Treatment	All Categories**
Max. user weight	200kg / 31.5st
Bariatric max. user weight	400kg / 63st
Medical Device Classification	Class IIa

*With clinical assessment. Support surface selection requires holistic patient assessment and cannot be decided on risk assessment score alone (EPUAP 2015).

**Both SDTI and Unstageable need to be treated as category 3-4 until depth of damage is established.

Mattress size

Series3 Auto (88 x 200 x 18cm)	PR13-760
Series3 Auto (105 x 200 x 18cm)	PR13-067
Series3 Bariatric (120 x 203 x 26cm)	PR13-761

Order Code