

OSKA®

Series3 Advance

Auto-weigh alternating air mattress



Alternating air, full replacement pressure care mattress.




Auto weigh



Quick inflation



Rapid repair feature

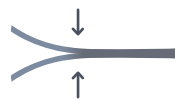

Welded seams



Zoned Offloading

Series3 Advance



Auto-weigh, active therapy alternating air mattress with quick 'release and replace' mattress covers.

The Series3 Advance is a full air mattress replacement system suitable for use across NHS Trusts.

The Series3 Advance has been designed with both clinicians and decontamination teams in mind. Its integrated air cell cover system reduces the level of decontamination required, enabling rapid turnaround between patients. Both the top and bottom covers can be quickly removed and replaced with clean or new ones, without disconnecting any air cells or tubing. This allows the mattress to be returned to service immediately, strengthening infection control protocols while enhancing operational efficiency in the acute care environment.

Effective pressure ulcer prevention & management

The alternating 1-in-2 cell cycle achieves regular periods of pressure redistribution and tissue offloading across the patient/support surface interface.

The auto-weigh pump offers a range of therapy settings and automatically adjusts to the patient's weight to ensure optimum air cell pressure at all times.

Bespoke zero pressure zones can also be created by using the unique 'quick connectors' to deflate specific air cells.



Infection prevention

Infection prevention is optimised by combining the following features;

- **Advanced Cover Technology** cover coating with inherent antimicrobial properties. Enhances physical robustness, doesn't crack and break down over time and reduces the impact of wear and tear on the cover, thereby minimising the risk of fluid strike-through into the underlying air cells.
- **ToughWeld seams.** Minimise the risk of fluid ingress into the underlying air cells.
- **Easyclean covers.** Can be machine washed at 95°C and tumble dried. They are also resistant to 10,000ppm available chlorine.



Enhanced patient safety & operational efficiency

The permanently inflated lower portion of each air cell reduces the risk of patients bottoming out during a power cut.

'Quick connectors' facilitate isolation and removal of damaged air cells without interrupting pressure area care therapy. This eliminates 'down-time' due to air cell damage and enhances efficiency by enabling individual air cell replacement in under a minute.



Series3 Advance



Giving mattresses a longer life!



Pressure relief



Patient comfort



Antimicrobial finish



Wipe clean



Machine washable
up to 95°C



Bleach cleanable
10,000ppm



Durability



Waterproof



Rapid
decontamination

Specification

Prevention Very High Risk*

Treatment All Categories**

Max. user weight 250kg / 39.3st

Therapy Modes Max Firm / Alternating /
Continuous Low Pressure /
Alternating Low Pressure /
Seat Inflation

*With clinical assessment. Support surface selection requires holistic patient assessment and cannot be decided on risk assessment score alone (EPUAP 2015).

**Both SDTI and Unstageable need to be treated as category 3-4 until depth of damage is established.

Mattress size

87 x 200 x 20cm

Order Code

PR13-420

Available for Rental - Hire - Outright Purchase